

Sterisafe DURO A3 Universal Steam

The intelligent system
for steam sterilization

Sterisafe DURO A3 Universal Steam

More efficiency and effectiveness for steam sterilization



Safe and sterile


Ideal for the theatre and other clinical applications

- Sterisafe DURO A3 Universal Steam is an innovative sterile barrier and packaging system in conformity with DIN EN ISO 11607
- Ideal for steam sterilization of rigid instruments
- Safe storage for rigid endoscopes
- The set comprises the Sterisafe DURO A3 Universal Steam and the Toolsafe Universal Lattice Base

Less is more

Less expenditure, more user-friendly

- No external wrapping required
- Lower costs for disposable and materials
- No seal to change because of the proven Pasteur's loop
- No maintenance required
- Reduced weight, easy handling

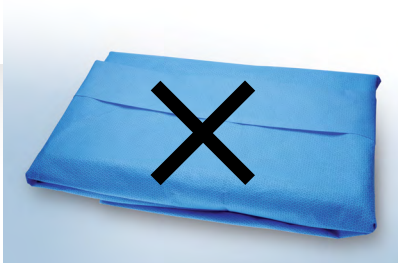

No time-consuming repackaging required.

Mehr Sicherheit, weniger Umstand


- High level of process certainty
- Suitable for steam sterilization of rigid endoscopes and instruments
- Integrated process indicators on the filter material
- Safe storage and transport
- For reuse the Sterisafe DURO A3 Universal Steam can be reprocessed in alkaline detergents at $> \text{pH } 10$
- Stockable

All applications under control

The Toolsafe Lattice Base A3 Universal made of silicone permits safe storage of different instruments.

Intelligent seal principle

Sterisafe DURO A3 Universal Steam operates on the principle of **Pasteur's loop** between the lid and the base. An additional silicone seal is no longer necessary.

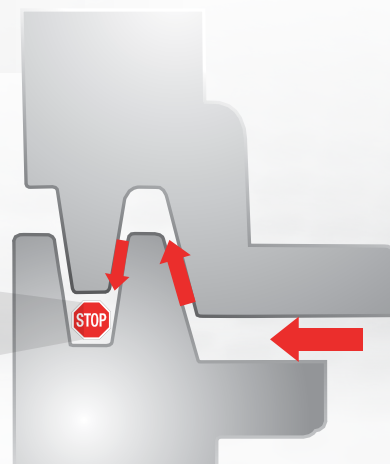


Hard-wearing

The use of translucent, tough, high-performance plastic makes Sterisafe DURO A3 Universal Steam extremely robust.

Sterisafe DURO A3 Universal Steam – so that everything runs smoothly

Rigid endoscopes and instruments are taken care of in every process phase. You can rely on that.

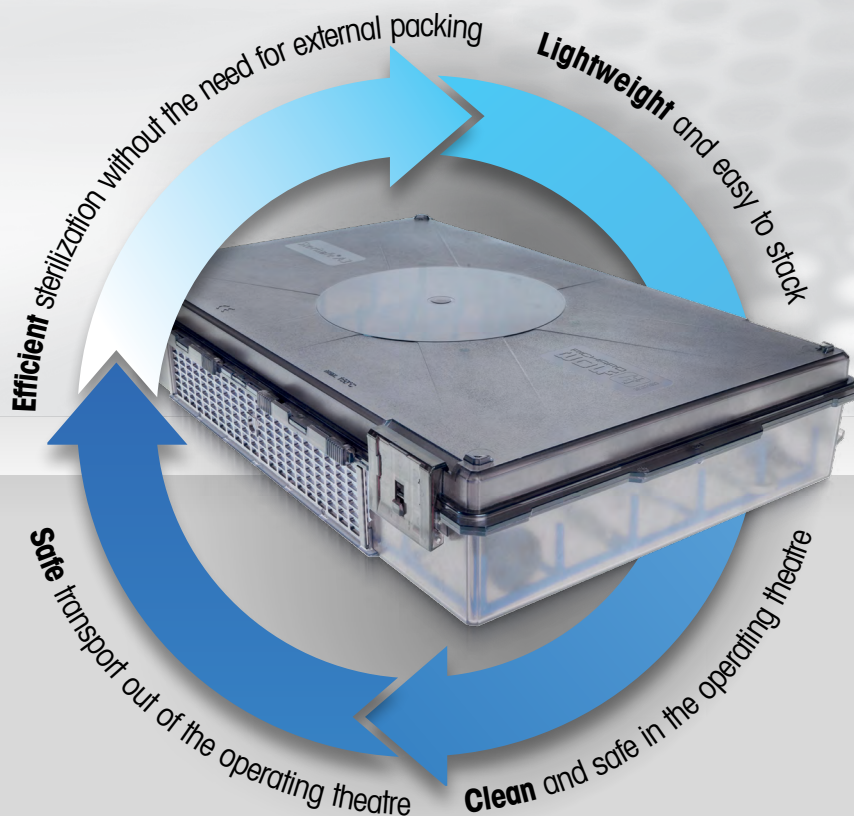

Sterisafe DURO A3 Universal Steam

and accessoires




spirit of excellence

Sterisafe DURO A3 Set Universal Steam, for steam sterilization

comprising:

Sterisafe DURO A3 Universal Steam
(382031003), Toolsafe Lattice Base A3
Universal (382031102), One-way filter
for steam sterilization (pack of 100,
382031501) 382032300

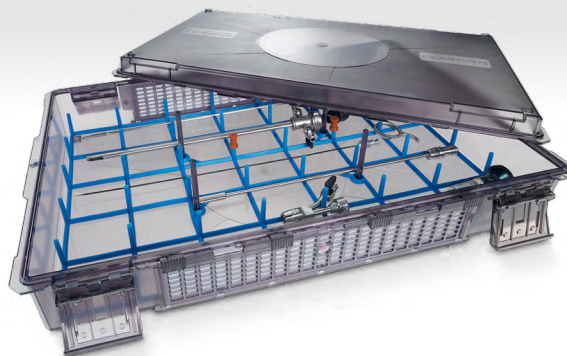

Accessories:

Single use filter with process indicator,

pack of 100,

for steam sterilization 382031501

Single use security seal,

pack of 100 382031202

Sterisafe Universal label holders,

universal 382031203



**For clear identification of the unopened status
of Sterisafe DURO A3 Universal Steam. The
security seal indicator is destroyed on opening.**



**For accommodating standard commercial
labels.**

Technical information:

Outer dimensions (L x W x H): 450 x 292 x 87 mm

Inner dimensions (L x W x H): 420 x 265 x 75 mm

Filter system: Single-use filter for steam sterilization procedure